

---

## Product Specifications

### LBL-Dye A832 NHS Ester

<https://proimaging.fr/home/41-lbl-dye-a832-nhs-ester-0210706.html>

### Ordering Information

Product Number: 0210706

Product Name: LBL-Dye A832 NHS Ester

Product Reference: CJ431

Unit Size: 1mg / 5 mg

Storage Conditions: Freeze ( $\leq -15$  °C), Desiccated, Avoid Light

Expiration Date: 12 months upon receiving

### Chemical and Spectral Properties

Molecular Weight:  $\sim 1352$


Solubility:

-Isotonic Dextrose Solution

-DMSO

Max Absorption: 825 nm

Max Emission: 832 nm

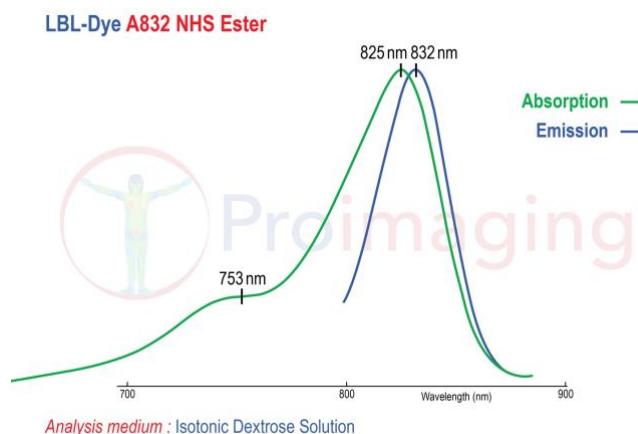

### Application Notes

LBL-Dye A832 NHS Ester is an activated bright dye with high photostability.

This dye is water-soluble and ideal for labeling sensitive proteins, peptides, and antibodies.

Product for Research Use Only. Not for Human Use.